Smart **Wi-Fi** body scale

User guide

Thank you for purchasing the **Smart Wi-Fi body scale** from Nexxt Solutions. We are excited to be part of making your life easier, safer and bringing peace of mind inside and outside your home.

Carefully read this manual prior to setting up and using this product.

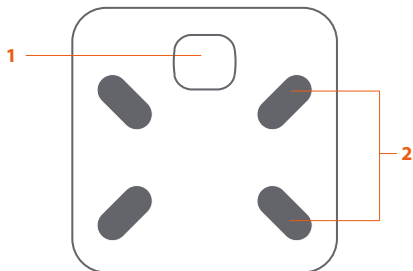
Package content

- Smart Wi-Fi body scale (1)
- Quick start guide (1)

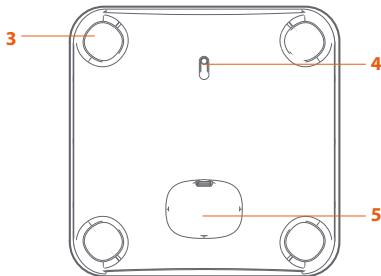
Product overview

1. LCD display
2. Electrode (each side, top and bottom)
3. Anti-slip feet
4. Unit reset button
5. Battery compartment

Top view

Bottom view

LCD display indicators

Symbol	Description
	Low battery indicator
%	Percentage

Symbol	Description
BMI	Body Mass Index
ft	Feet
cm	Centimeter
kg	Kilogram
St	Stone
lb	Pound
	Wi-Fi connection
	Data is ready to be transmitted to the app: <ul style="list-style-type: none">- If the data transmission fails, the symbol will remain. The pending measurement data will be temporarily kept in the scale and transmitted to your smartphone when the next measurement is complete.- The symbol will disappear after successful data transmission, and you may check your personal health data stored in your smartphone.

Preliminary steps

Battery installation

1. Turn the scale over and place it on a hard surface.
2. Open the battery compartment on the bottom of the scale.
3. Insert three AAA batteries (not supplied) observing the correct polarity inside the compartment.
4. Replace the battery cover and wait until **0.0kg** are displayed on the screen.

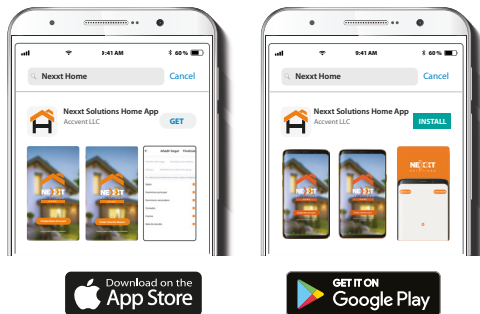
When the “” appears on the screen, the device will power off in about 4 seconds. That means the batteries are low.

All three batteries must be replaced at the same time.


Do NOT mix new batteries with the old ones.

Installing the Nexxt Solutions Home App


Download Nexxt Home app from Google Play or the Apple App Store.

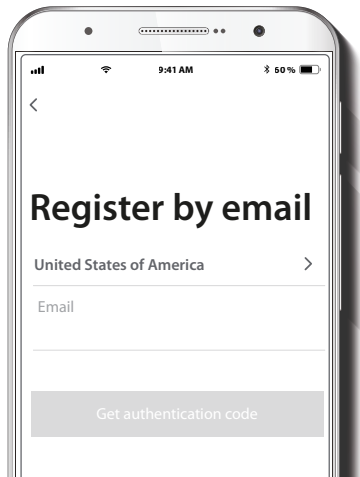Step 1


Register by email. Fill in the required fields.

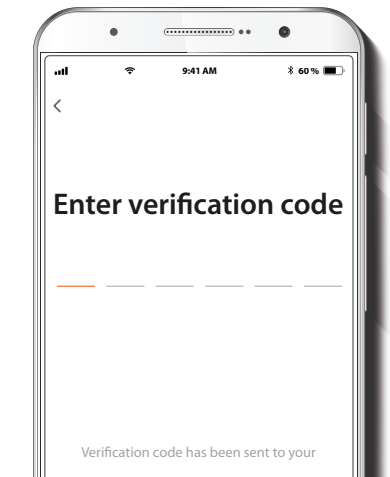Step 2


Enter the verification code you received.

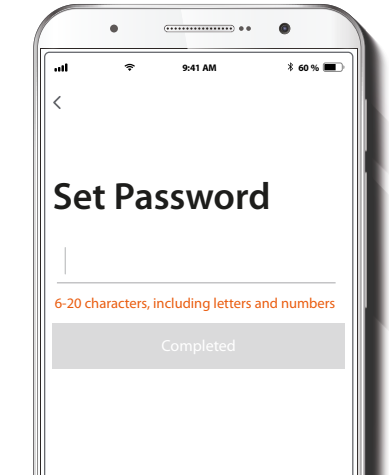Step 3

Create a new password and then select **Completed**.

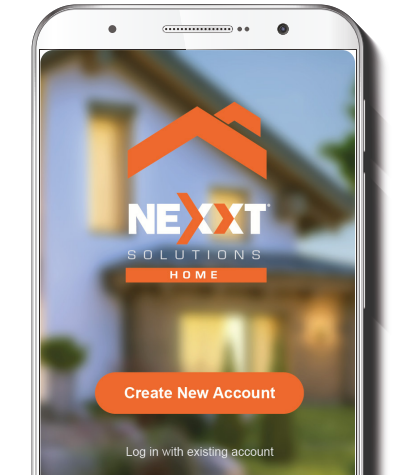
Step 4

The Nexxt Solutions Home account has been successfully set up.


Add device

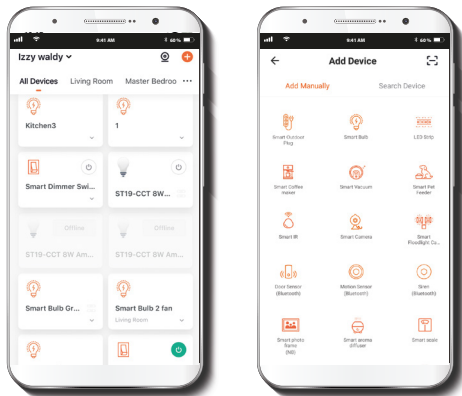Step 1

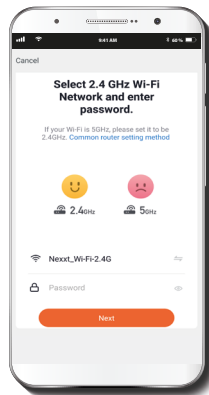Select the + icon to add the device, followed by the **smart scale**.


Step 2

You must connect to a **2.4GHz WiFi network**. Enter your network password. Make sure you enter the correct password, otherwise the device will not connect.


- Press and hold the reset button on the back of the scale.

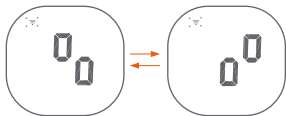• The  icon and the  icon will appear alternatively on the LCD, accompanied by a blinking Wi-Fi icon. This indicates that the device has entered into the pairing mode.


Note: If the Wi-Fi icon blinks slowly on the LCD, you can click the **Other mode option** on the top right corner to manually select **AP Mode** and pair the scale.

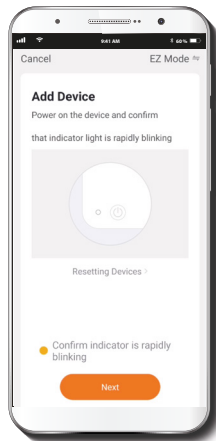Step 3

On the app, press **Confirm indicator is rapidly blinking**.

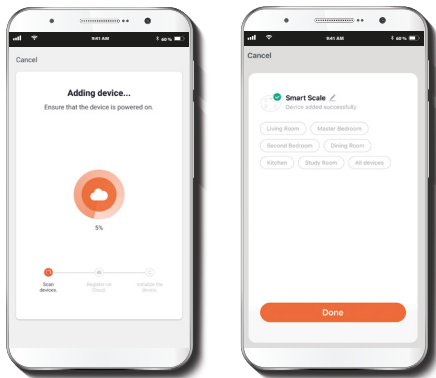Step 4

Device added successfully will be displayed. Click on **Done**. Once the pairing is complete, the LED will turn off.

Once the pairing is complete, the  icon will appear, making the Wi-Fi symbol become steady on the LCD.

If the pairing process fails, then the **Conf** notification will be shown instead.

Selecting the measuring unit

After adding the device successfully to your network, you can select the measurement unit by pressing the **unit reset** button located on the bottom of the scale. The available options are: kilogram, pound or stone as your measurement unit. This can be selected through the app as well.

Starting a measurement

Step 1

Step on the platform barefoot.

Step 2

Stand still and keep full contact with the electrodes until the LCD stops displaying the **0000** sequence.

Step 3

Make sure your Wi-Fi connection is active, so the analysis of your body composition can be sent wirelessly to the app.

Make sure to remain on the scale until the **Data transmission** symbol disappears from the screen.

Step 4

First tap on "View" message to claim your weight information. Then, press the **body weight results** to input your personal information using the app, and click latest measurement record to distribute your User ID.

1. Press the **body weight results** tab on the app (picture #1) to review your personal weight-in results for your body metrics (picture #2).

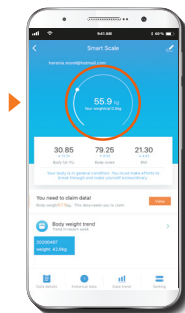
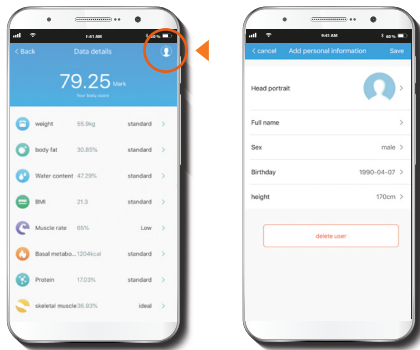
Imagen #1


Imagen #2

2. Click on the  icon to update your personal information.


Step 5

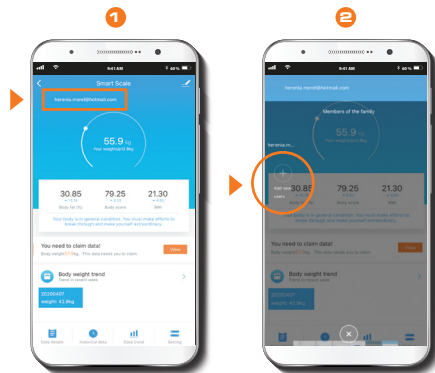Adding additional users:

1. On the app, press the upper left-hand side of the mobile screen. (By acct. #)

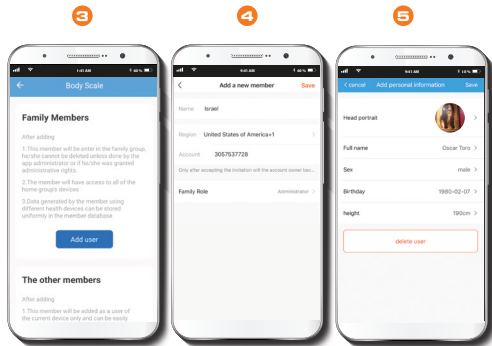2. Press the (+) symbol to **add a new user**.


3. Select **Member of the family** or **The other members** option, followed by the **Immediately add** button.

4. Enter the member's information and set under family role the permission granted to that user.


5. To delete a user, press and hold the name you wish to eliminate. Press **Delete** user to permanently remove it from the list.

Daily measurement

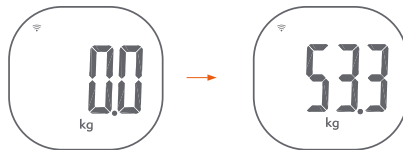- Stand still and keep full contact with the electrodes until the LCD displays the weight.
- Make sure your Wi-Fi connection is active, so the analysis results will be automatically matched with the user ID you created and send them to your phone through the app.

Application main menu

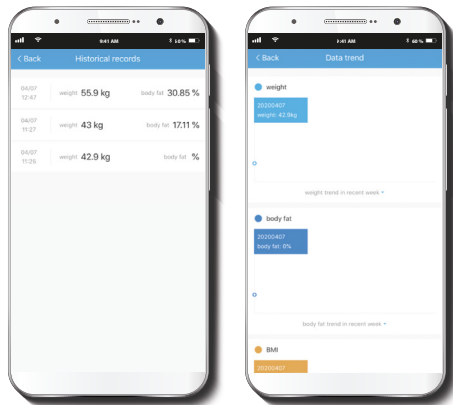
Historical data: allows you to maintain and track all measurement records on a weekly, monthly or yearly basis.

Data trend: provides a graphical representation of your weight, body fat and BMI (body mass index) on a weekly, monthly or yearly basis.

Setting: Allows you to select the unit of measurement you wish to use.

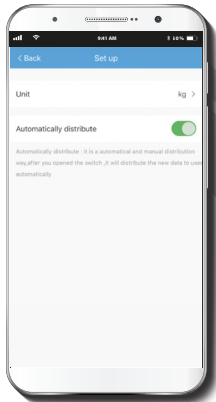Body Metric Data: provides an itemized report of all body related measurements.


Set up:

- **Unit:** allows you to set your preferred unit of measurement.
- **Automatically distribute:** use it to enable or disable the distribution of data to users automatically.

Settings

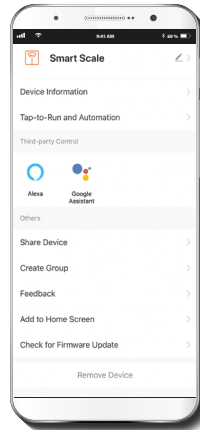Click on the pencil icon on the upper right-hand corner.

Device name: allow you to change the name of your device, icon image and assigns location

Device information: will show the information related to the Wi-Fi smart scale.

Tap-To-Run and Automation: links the Wi-Fi body scale with other smart products to automatically control them through a customized command.

Third-party control: enables the voice recognition feature through compatible devices.

Share devices: allows sharing the Wi-Fi body scale with family and friends.

Create Group: makes management of devices easier by selecting and arranging multiple components under a single group.

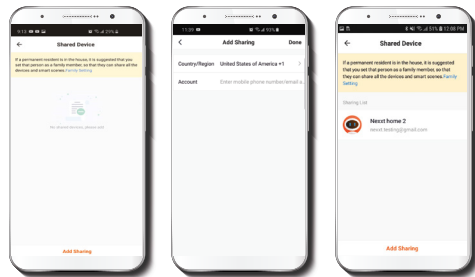
Feedback: let us know if you have any issues setting up the product.

Add to home screen: adds a shortcut of the device to the main screen on your mobile phone.

Check for firmware upgrade: we will sometimes add new features or provide updates with a new firmware upgrade.

Remove device: allows you to delete the Wi-Fi body scale from the application.


Share device

Click on **Share devices** then tap **Add sharing**. Enter the mobile number/e-mail address of family and friends. Tap on **Done**. The newly configured users will now have access to your smart scale.

Note: family or friends will need their own Nexst Solutions Home Account to access the device.

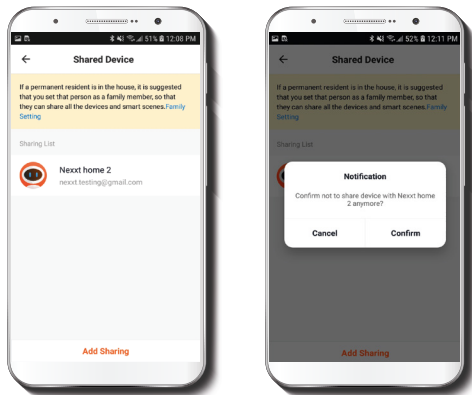Delete user

To remove a user, on Android press and hold the name you want to eliminate. On iOS, slide to the left the user you want to remove. Then click on **Delete**. A notification confirming you don't want to share the device with this user will be displayed. Click on **Confirm** to finish.

Troubleshooting

Error display	Description	Solution
	Overload. The device will power off	Stop using this scale for measurement
	Low battery. The device will power off in four seconds	Replace all three batteries at the same time
	Wi-Fi connection failure	Rejoin the wireless network

When measuring....

Issue	Root cause	Solution
Abnormal measurement results: - Too high - Too low - Huge difference between two recent measurements	Incorrect posture	Step on the platform barefoot and stand still
	The device is placed on soft carpeting or uneven surface	Place the device on a flat, hard surface
	Cold body that may result in bad blood circulation	Warm up your hands and feet to resume blood circulation and then measure again

Issue	Root cause	Solution
Abnormal measurement results: - Too high - Too low - Huge difference between two recent measurements	Cold electrodes	Place the device in a warm room for a while and then measure again
	Your feet are too dry	Wipe your feet with a damp cloth, keeping them slightly damp
Scale does not proceed to analyze body composition	Taking a measurement while wearing socks or shoes	Wipe your feet with a damp cloth, keeping them slightly moist
	The user fails to select the User ID from what the system found	Please assign a user ID

FCC ID#: X4YHASS130

Nexxt Solutions® is a registered trademark. All rights reserved. All other brand names and trademarks are the property of their respective owners. Made in China.

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